

# KONQUEST **KDT-2106v2**




## QUICK USER GUIDE

### **KONQUEST KDT-2106v2**


#### DIGITAL INFRARED THERMOMETER

# PRODUCT STRUCTURE

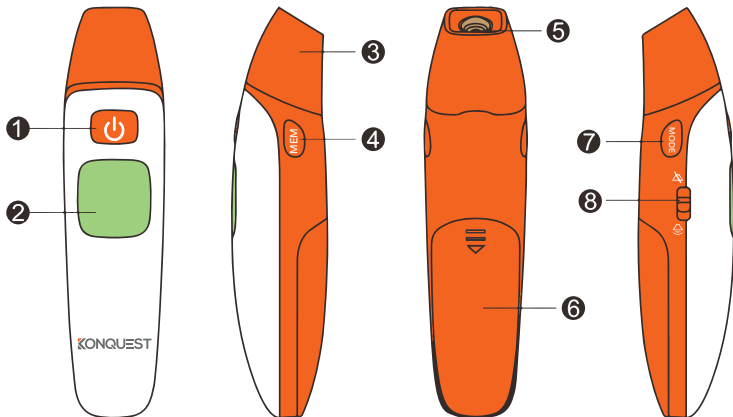


- [1] Power button / Measure button
- [2] LCD display screen
- [3] Probe cover (keep it when measuring forehead temperature)
- [4] Memory button / Unit switch button
- [5] Probe cover (remove it when measuring ear temperature)
- [6] Battery cover
- [7] Mode button
- [8] Sound switch



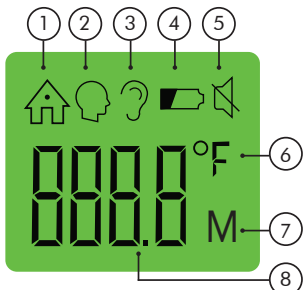

# DISPLAY DESCRIPTION

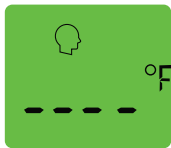


- [1] Object temperature mode
- [2] Forehead temperature mode
- [3] Ear temperature mode
- [4] Low battery
- [5] Mute / Un-mute
- [6] Fahrenheit / Celsius degrees
- [7] Memory recall
- [8] Temperature value





## 1 FOREHEAD MODE TEMPERATURE MEASURE

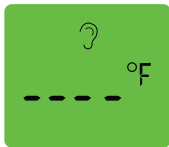


- [1] Put the cover on the probe.
- [2] Press and release the **POWER/MEASURE** [  ] button for one second.
- [3] Press the **MODE** [  ] button. Select **FOREHEAD**. The symbol **HEAD** [  ] is displayed in the screen.
- [4] Place the probe with the cover on the forehead between the eyebrows.
- [5] Press the **POWER/MEASURE** [  ] button to start measurement.
- [6] After the beep temperature will appear in the LCD screen.

## 2 EAR MODE TEMPERATURE MEASURE



- [1] Take off probe cover.
- [2] Press and release the **POWER/MEASURE** [  ] button for one second.
- [3] Press the **MODE** [  ] button and select **EAR**. The symbol **EAR** [  ] is displayed in the screen.
- [4] Insert the probe in the ear canal.
- [5] Press the **POWER/MEASURE** [  ] button to start measurement.
- [6] After the beep temperature will appear in the LCD screen.

### 3 OBJECT MODE TEMPERATURE MEASURE

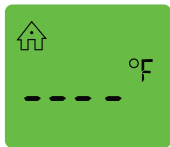

[1] Put the cover on the probe.

[2] Press and release the **POWER/MEASURE** [  ] button for one second.

[3] Press the **MODE** [  ] button. Select **OBJECT**.

The symbol **HOUSE** [  ] is displayed in the screen.




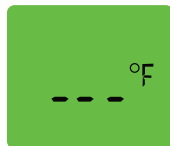
[4] Point the thermometer to the center of the forehead between the eyebrows.

[5] Press the **POWER/MEASURE** [  ] button to start measurement.

[6] After the beep temperature will appear in the LCD screen.

### 4 SWITCH FAHRENHEIT TO CELSIUS DEGREES

[1] Press and release the **POWER/MEASURE** [  ] button for one second.



[2] Press and hold the **MEMORY/UNIT** [  ] button for more than 4 seconds. °F or °C are flashing.

[3] Press the **MEMORY/UNIT** [  ] button to switch between °F or °C.

[4] The thermometer shuts down after 10 second.  
The setup is successful.

### 5 MUTING THE THERMOMETER



[1] Toggle the **SOUND SWITCH** [  ] to switch between MUTE and UN-MUTE.

[2] The [  ] symbol is displayed in MUTE MODE and disappears in UN-MUTE MODE.

READ INSTRUCTION MANUAL THOROUGHLY BEFORE USE